
Nikon N-STORM Super Resolution Microscope





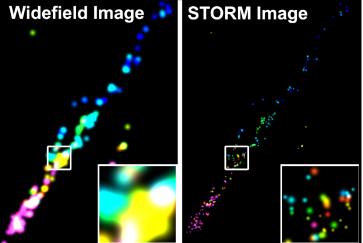

Instrument Description

The Nikon "Storm" system is an N-STORM super resolution microscope, capable of achieving 20nm XY and 50nm Z resolution - 10X increase in resolution over a traditional fluorescence microscope. STochastic Optical Reconstruction Microscopy (STORM) is our most advanced imaging modality offered at the MSR and is ideal for experimental questions requiring nanoscale resolution. This instrument is capable of acquiring up to 3 fluorescence channels and has the optics to acquire 3D STORM images. Due to the advanced nature of STORM imaging, the MSR offers an Assisted STORM Imaging Service to assist users with no background in STORM. This service includes the preparation of MSR supplied STORM buffers and MSR staff performing the STORM imaging with user guidance.

Instrument Features

- Super resolution imaging (20nm XY and 50nm Z resolution)
- Optics to acquire up to 3 fluorescence channels (488nm / 561nm / 640nm)
- 3D STORM imaging (with astigmatism lens)
- Photoactivated Localization Microscopy (PALM) imaging for samples with photoconvertible fluorescent proteins

Nikon N-STORM Super Resolution Microscope



Instrument Details	
Available Lasers	405nm / 488nm / 561nm at 100mW 640 nm at 125mW
Available Objective Lenses	20x Plan Apo (NA = 0.75; WD = 1mm; DIC) 100x Oil SR Plan Apo (NA = 1.49; WD = 120μm; DIC)
Available Camera	Andor iXon Ultra 897 - Back-illuminated EMCCD sensor - Pixel size = 16µm - Sensor Size = 512 x 512
Emission Filters	DAPI - 455nm / 50nm GFP - 525nm / 36nm Cy3 - 605nm / 52nm Cy5 - 705nm / 72nm
Supported Samples	IMPORTANT: Samples must in the format below Glass-bottom Dishes (35mm) with No. 1.5 coverslips Recommended No. 1.5H coverslips Mattek P35G-0.170-14-C Unmounted samples only (store in 1x PBS)